
SIL PROMESH® **SIL PROMESH**

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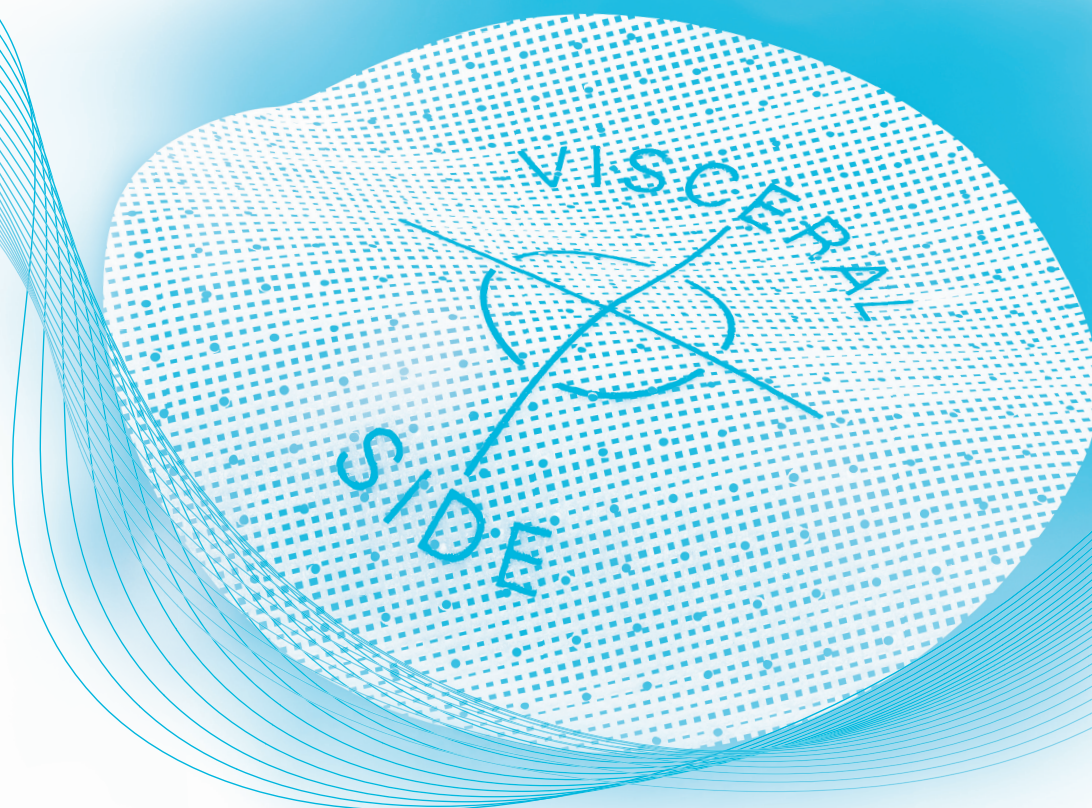
The optimal combination of polypropylene and silicone


- One side made of non woven polypropylene (50 g/m²) to allow an optimal tissue integration
 - Development of a fibrous anchoring by the cells of the peritoneum in the non woven polypropylene matrix.
 - Pore size designed for a good tissue integration, and a limited risk of seroma formation^[1].
- One side with a visual marker made of non adhesive silicone to avoid the formation of visceral adhesions
No tissular colonisation on the hydrophobic silicone film.

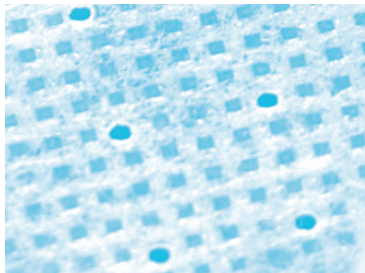Struktura matrycy polipropylenowej
Structure of the non woven polypropylene matrix

Perforowana, antyadhezyjna warstwa silikonu
Non adhesive silicone film with macroperforations

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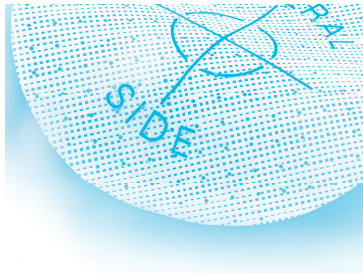Optimal physical characteristics

- Two intimately bounded sides:
Specific assembling process to avoid any risk of delaminating, even after cutting the prosthesis.
- Macroperforations:
Seroma drainage of the surrounding tissues.

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Widok na otwarto
Open surgery view

Safe and easy implantation

- **Excellent shape memory and flexibility:**
Easy deployment and easy handling by laparoscopy or open surgery.
- **Safe covering of the defect:**
Target mark on the silicone side guiding the positioning of the translucent prosthesis.
- **Ready to use, without prior preparation.**
- **Compatible with usual fixation devices:**
Sutures, staples, helicoidal fixations, glue... absorbable or not.

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Long term clinical efficacy and tolerance Prospective multicenter clinical evaluation

- 115 patients prospectively included and followed-up
- 36 months follow-up
- 7 French centers
- Detection of adhesions between the prosthesis and the viscera by ultrasound examination^[2] performed since the 6th postoperative month in 60 patients
- Intermediate results published^[3]

A proven anatomical efficacy:

- Low recurrence rate
- No recurrence incriminating the prosthesis

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
A range adapted to all indications
and surgical approaches

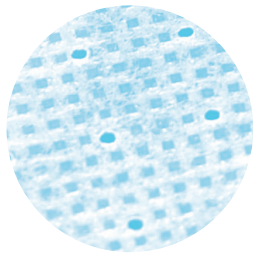Kod produktu Reference	Rozmiar [cm] Size [cm]	Kształt Shape	Zastosowanie w laparoskopii Laparoscopic use
SI09RMP	Ø 9		X
SI12RMP	Ø 12		X
SI18RMP	Ø 18		X
SI1015MP	10 x 15		X
SI1515MP	15 x 15		X
SI1520MP	15 x 20		X
SI2020MP	20 x 20		
SI2030MP	20 x 30		
SI2530MP	25 x 30		
SI3030MP	30 x 30		
SI3846MP	38 x 46		

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